


OPERATING INSTRUCTION

Household Beauty Instrument
Ultrasonic Ion Shovel Skin Instrument

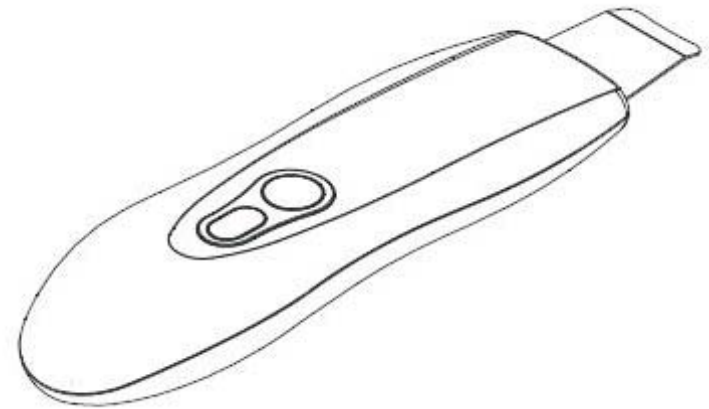

Thanks for choosing our product.
Please read the operating instruction of this instrument
for correct and safe use, so that reach to the best effect.
If this product is lent to others, please enclose the
operating instruction and keep it properly.

CONTENT

Product Name	2
Product	2
Accessories	2
Specification Parameter	3
Charging Method	3
Three Modes	5
Using Method	6
Primary Key Points	6
Confirmation for Ultrasonic Anion	6
Normal Use	7
Skin Care Skills	9
Key Points of Skin Care	9
Face Care	9
Cheek Care	10
Nose Care	10
Product Maintenance	11
Problem Solution	11
Trouble Shooting	12
After-Sales Service	12
Warranty Rules	12
Quality Warranty, Warranty Card	13

WARRANTY CARD

No.			
Item Name		Type	
Product Name			
Customer Name		Tel.	
Address			
Purchasing Date			
Free Repair	Object Part	Complete Instrument	Except for Consumables
	Time	1 Year	From the Date of Purchasing Date
Store Name			
Record of Repairs I			
Record of Repairs II			
Record of Repairs III			
Warranty Rule	See the Warranty Notes		

QUALITY WARRANTY

1. If you use it in normal way in accordance with the notices of operating instruction and label enclosed to beauty instrument, for the failure, you go to the sale point for free repair.
2. During warranty period, if the failure happens, please go to the sale point for free repair with this instrument and warranty certificate.
3. Please tell the address in detail for return.
4. The gift is not in free repair range.
5. During warranty period, for the following situations, need pay for repair.
 - A. The failure and damage due to wrong method, wrong repair and refit under your own disassembly.
 - B. The failure and damage due to dropping on the floor after buying.
 - C. The failure and damage due to fire, earthquake, flood and other natural disasters.
 - D. The failure and damage of lithium battery, consumable and other parts.
 - E. The failure and damage due to transport, carrying and dropping after buying.
 - F. The situations that are not listed on warranty certificate.
 - G. The situations that there are no date, customer name and sale point name or seal on the warranty certificate.
 - H. This quality warranty is only effective within China domestic market.
6. The free repair is only for the beauty instrument which has clear time and meets the condition. For other situations, it shall have certain restriction. For the repairs after warranty expires, please contact with the seller.
(Such as the replacement of battery and consumables.)

STATEMENTS:


For the better service, this operating instruction may be changed according to the update of product, do not notice.

If the actual operating of the product does not comply with the instruction, please access to the website of this company for information of new product or contact with the customer service. This company is not legally responsible for the losses due to this instrument or various severe requirements, and this company reserves the right to interpret all above.


PRODUCT NAME

Product

Product Appearance & the Name of Each Part

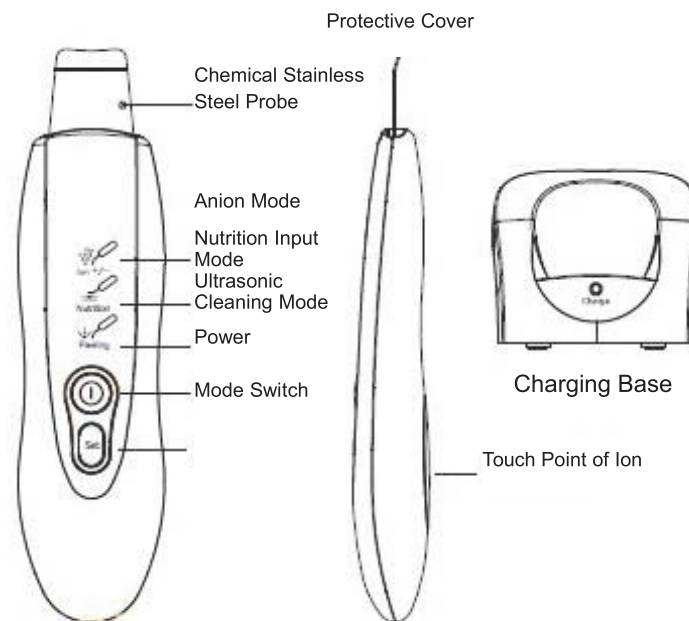

Accessories

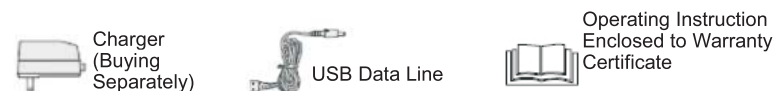

This product supports environmental protection concept, it can be charged by common mobile phone charger (5V, ≥500mA). Accessories only include USB data line without charger, please buy it separately if the user need charger.

Product Parameter

Specification Parameter:

Battery	Product Size /mm	Produce Weight
Lithium Battery	About 180*45*22	About 90g

Wireless Induction Charging:

Charging Voltage	Charging Current	Charging Time	Using Time
DC 5V	About 200mA	About 3-4H	>1H

※The using time is different if different loads are used in daily.

Charging Method

Notice the charging order, connect to the power supply firstly and then put the product!

- 1) The charger is inserted to AC socket, the USB terminal is inserted to charger, the DC terminal is inserted to charging base, if the indicator light of charging base is on, that indicates it is in normal working.
- 2) Turn off the product to put into the product base vertically, keeping the close contact of bottom, that shall produce best charging effect.



Trouble Shooting

If there is abnormal during product using, please confirm whether there are the following phenomena:

There is no reaction when pressing the power key.

It is estimated that its power off, please charge it.

No spray, or less spray

The spray is difficult to emit, if there are water or greasy dirt around shovel blade center. Please wipe shovel blade with clean cloth.

Charging lamp is off.

The plug on charging bench is not connected well. Please connect it again.

Sometimes due to the full power, after use it for 10 minutes, and then charge it and confirm the lamp state.

It can only be used for a while after charging.

Confirm whether the charging bench is placed well.

After-Sales Service

Charge battery is consumable (It is not warranted).

It is only used for short time after charging, which proves the battery comes to its service life (Normal charging times are 300 times).

For the replacement of charge battery, please contact with store or after-sales service point.

Warranty Rules

For the failure during normal using, one year warranty (except for consumables). Within warranty period, any repair fee is not needed but transport charge. For the failure, please contact with seller or after-sales service point.

SKIN CARE SKILLS

1. Press the Power key, and turn off the power.
2. Use clean cloth to wipe out the drop of water on machine.
3. The front mental part of machine shall be cleaned by dry cloth, but when it is very dirty, use cloth to soak neutral cleaner, then use dry cloth to clean it.

ATTENTION: Forbid to add chemical agent, like cleaning agent as thinner, volatile oil and grinding agent.

4. The machine will be functionless without cleaning the dirty. So it must be cleaned after using.
 5. After using, zinc plate must be covered with protection cap, in case of robbing off drop of water.
- ※ **The machine itself has handled with drip-proof, but charger, charging base etc, need water-proof.**
- ※ **The process of drip-proof aims at spray, different from water-proof. So it cannot soak in water or wash.**

PROBLEM SOLUTION

Problem Solution :

1. Do every time three modes need be all used? How to decide?
The three modes have different functional modes. Choosing one mode on basis of one need is all right.
2. Can it be used on other place apart from face?
The product is developed for face, nose care etc. Suggesting not use to other place.
3. How many times can it be used one day?
The normal care is a time per day, but 2-3 times a week is better. Do not use it too often.
4. Why the skin care place must be moistened?
Peeling mode and Nutrition mode aims to spread supersonic wave to skin. The beauty effect of shovel blade in dry part is less than water part. Suggesting not use without water.
5. Can it be used in dry condition?
Yes, it can. Please use LON-anion mode. After cleaning face, anion is very suitable.
6. Can soft skin product replace water?
Water soluble product can be used, but avoid using product which strongly stimulate skin. Especially for sensitive skin, suggest use our recommend product.

Charging Method



It is not allowed to charge at bath room or moist site. Although the product has drip-proof function, the plug and charging base do not have.

- 3) Charging need last for about 3-4H, 4H are required for the first use. In daily, it is ok to charge fully. With the charging is completed, the lamp on the product shall be off, but the lamp on base shall not be off. After the charging is completed, the plug shall be pulled out from power supply, as long electrifying time shall waste electricity! And the service life of charging base shall be affected. If the using method of product is different, the service life of product is also different.



With power, it is prohibited to place metal or other conductors on the base. It is prohibited to connect to the power for long time, and the power plug shall be pulled out after use.

- 4) Please turn off the power, if the power lamp blinks during use. It is commended to use special charging base to charge. If the noise appears on TV and radio during charging, additional external socket shall be used. Charging under lower temperature may reduce the battery service life. Please charge under appropriate temperature. The appropriate temperature is between 15-30°C. It is not allowed to charge under sub-zero temperature. A little heat emitting during charging is normal rather than instrument failure. In order to protect the built-in battery, please charge for 1 time per month, 3-4H per time, if the instrument is left unused for long time. (If the power is full, need not charge again.)

THREE MODES

Three Modes

Three modes are three different functional modes! One mode shall be selected as your demand.

A. Peeling Mode

If select this mode, the corresponding Peeling green lamp shall be on. The main effect of this mode is deep clean. The greasy dirt shall be eliminated under the ultrasonic effect. If the greasy dirt residual appears on metal part, please clean it with clean towel or tissue.

B. Nutrition Mode


The main effect of this mode is to input the nutrition, tap and massage. Shovel shall be used by flatting on the skin under this mode, which can tap surface and help to nutrition input. If select this mode, the corresponding Nutrition green lamp shall be on.

C. LON- Mode

If select this mode, the corresponding LON green lamp shall be on. The shovel can be used by flatting on the skin or in 30 degrees. Its using direction is opposite to that of Peeling mode, pulling from up to down. The contacting area shall be adjusted according to your need! It can realize ion skin rejuvenation and phytosome.

Cheek Care

- 1) Chin to ear.
- 2) Mouth to temples.
- 3) Nasal side to temples, chin, mouth corner, skin care of cheek surrounding.

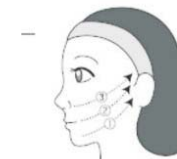

- 4) Chin to mouth corner. (Ch Mouth Corner)
- 5) Mouth corner to zygoma. (Cheek Care)



- 6) Chin to bridge of nose. (Mouth Corner Care)




ATTENTION: Please use the low-grade for angulus oris care.

The above steps from 1 to 6 are combinations, please crossover right to left for five to ten minutes.

※The machine will automatically stop in the first ten minutes.

Nose Care

- 1) Wing of nose, following the arrow from bottom to top slide the zinc plate on nasal side.
- 2) Slide to the front of arrow.
- 3) Slide from bottom to top on sides again.
- 4) Slide slowly from bottom to top.

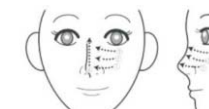

1 to 4 of skin care is a whole and needs 2 to 3 minutes.(If nose pore is large, one shall add time)

- ❗ 1) Please do not slide on eyeball, eyelid, ear and neck.
2) Please do not slide on the place where eczema, inflammation, acne and pimple can occur. When one is in or after skin care, the skin is uncomfortable like redness, measles and itch, one shall stop use it. And use towel for cold compress. Then smear the daily skin care product.

SKIN CARE SKILLS

Five key points of skin care

When the skin is warm, it is more effective for skin care. After taking shower, using warm towel to apply face or using warm and cold beauty instrument, etc., the temperature of body is warm.

Spray employs clean water turning into mist. Toward the place of care to spray, then starting the machine to use but the zinc plate cannot stand too much strength. (If your skin is dry, there is no respond; so please add water in time.)

Please move the zinc plate at the speed of three second of ten cm. (Avoid pausing at one place for a long time)

The time of skin care is one time per day (about 5-10 minutes), but 1-2 times per week is better.

After using the machine, the method of skin care is smearing toning lotion, latex etc. Thus to add water.

Cleaning skills

The zinc plate need to slide slowly. (not shaving forcibly, but sliding gently)

Greasy dirt will show up by vibrating.

One time drops enough water for the same place. (If you have skin allergy, please do not slide too much in the same place.)

※If greasy dirt can not show up for one time, regularly use it can wipe it out.

Face Care

1) Sliding in straight line from top to bottom.

2) Sliding toward the front along marker.

3) Sliding from bottom to top again.

※Not shaving forcibly, but sliding gently.

※If you have skin allergy, please do not slide too much in the same place.

※LON-mode and Peeling mode should use in a reverse direction.

The following will introduce care method of Peeling mode. The care method of Nutrition and Peeling is the same.

USING METHOD

Confirmation for ultrasonic effect:

Before use, please take off the protective jacket firstly, and then confirm it.

The skin is wetted with water or sprayer, use the shovel to touch the skin slightly and slide it on the skin slowly. If the spray appears at the front end of shovel, it indicates that ultrasonic occurs. When the spray is less, you can increase the water quantity. (Because of the light and background color, sometimes it is difficult to see the spray. So you can feel it with hand and rapid microvibration is normal)

Key points of primary test:

※ This instrument does not have auto power-off function, which shall stop automatically after running for 5 minutes. Product's temperature is a little higher, that is normal phenomenon.

※ If the auto power-off happens or the shovel is weak it indicates the voltage is not enough. Please turn off the power to charge.

※ Please confirm whether it is charged well before use. The protective jacket of shovel shall be taken off firstly, and then use it.

※ The skin care must be done after makeup removal.

※ It is best to use 1 time per day, 1-2 times per week. Over frequent use shall cause the discomfort of the skin.

※ When use it at the first time, please test your skin firstly. Clean the place outside of the face and the water is wiped by the towel.

Make a little water on the cleaned place, and start to use the instrument at the small area after connecting to the power. If there is no abnormal 48 h later, you can use it. If there is abnormal please stop using.

Confirmation for anion:

Contact the conducting position of rear cover with hand, and face flat on the shovel. If anion generates the acupuncture feeling, it works well.

Hand Touch Point of Anion



1. If there is acupuncture feeling on anion, for the new user, please hold the rear cover slightly to adapt to it slowly.

2. The new users or users with sensitive skin shall use less time. The state of the skin shall be observed while using. Check whether there are measles, pustule, getting red or itch etc., if any, please stop using.

NORMAL USE

※ Confirm whether it is charged well prior to use. Please take off the protective jacket of shovel firstly.

1) Please remove your makeup prior to use.

2) Take off the protective jacket.


3) Spray water on the face (Better to use the sprayer).

(Although the water-soluble hair gel can be used, massage cream has bigger stimulation on the skin. So try to avoid using it.)


4) Keep pressing the POWER key  to enter into standby mode.

The three LED lights Peeling / Nutrition / LON-  all blink.


5) Press the SET key  to select the mode:

6) Press SET key, select Peeling mode , the corresponding Peeling green light shall be on. Main effect is to clean deeply, discharge the greasy dirt under ultrasonic wave.



7) Press SET key, select Nutrition mode , corresponding Nutrition green light shall be on. Shovel shall be used by flatting on the skin under this mode, which can tap surface and help to nutrition input.



8) Press SET key, select LON-mode,  the corresponding LON-green light shall be on. Shovel shall be used in 30 degrees, and then it can be adjusted according to your demand. It can realize ion skin rejuvenation and phytosome.

9) Press SET key to back the standby mode.

10) After one cycle, you can keep pressing the POWER key to turn off. You also can continue to press SET key to select the mode you like to use continually.



KINDLY REMINDER:

LON is also called as electric face-lift, with acmesthesia, which is normal phenomenon. Please adapt to it slowly.